



User Manual

Utility Chair & Desk Utility Chair

SEATING SOLUTIONS



Intended Use

Thank you for choosing a **K Care Utility and Desk Utility Chairs**.

These chairs are intended for sitting in a homecare, aged care or hospital environment. They are height and width adjustable to suit users of various size and mobility. The Desk Utility design is intended to allow users to position themselves closer to desks and tables.

Specifications

KA586KDV03 Champagne Utility Chair

Seat Width	460 – 510 mm
Seat Depth	460 mm
Seat Height	430 – 580 mm
Overall Depth	510mm
Backrest Height	440mm
Max User Weight	200kg

Materials: Mild steel with powder coat finish and padded upholstery covered in vinyl.

KA586KDDV03 Champagne Desk Utility Chair

Seat Width	460 – 510 mm
Seat Depth	460 mm
Seat Height	430 – 580 mm
Overall Depth	510mm
Backrest Height	440mm
Max User Weight	200kg

Materials: Mild steel with powder coat finish and padded upholstery covered in vinyl.

KA586KDG05 Grey Stone Utility Chair

Seat Width	460 – 510 mm
Seat Depth	460 mm
Seat Height	430 – 580 mm
Overall Depth	510mm
Backrest Height	440mm
Max User Weight	200kg

Materials: Mild steel with powder coat finish and padded upholstery covered in vinyl.

KA586KDDG05 Grey Stone Desk Utility Chair

Seat Width	460 – 510 mm
Seat Depth	460 mm
Seat Height	430 – 580 mm
Overall Depth	510mm
Backrest Height	440mm
Max User Weight	200kg

Materials: Mild steel with powder coat finish and padded upholstery covered in vinyl.

Safety Precautions


Safe Use Environment

This product is safe for use in dry environments (i.e. not to be used in showers, baths or other wet area environments).

This product may be used outdoors in dry conditions for short periods however should not be left exposed to the elements for extended periods of time or degradation of the upholstery may occur.

Safety Warnings

- Do not exceed the maximum user weight
- Periodically inspect the assembly fasteners to ensure they are tight
- Ensure that the leg adjusting push pins are fully engaged before use:

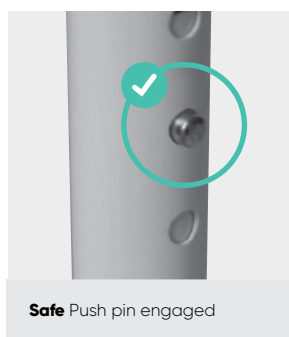


- Ensure that the chair is stable and level on the floor before sitting
- Ensure there is adequate space around the utility chair to enable safe movement onto, off from and around the chair

Assembly Overview

After unpacking the box you should find the parts listed below.
If there are any missing parts, please contact:

K Care Customer Service: 1300 783 783 / customerservice@kcare.com.au


- | | | | | |
|---|-----------------------------------|----|---|---------------------|
| 1 | Utility Chair Side Frame x 2 | OR | 5 | Spacers x 4 |
| 2 | Desk Utility Chair Side Frame x 2 | | 6 | Handwheel Knobs x 8 |
| 3 | Seat Frame | | 7 | Allen Key |
| 4 | Padded Backrest | | 8 | Screws x 4 |

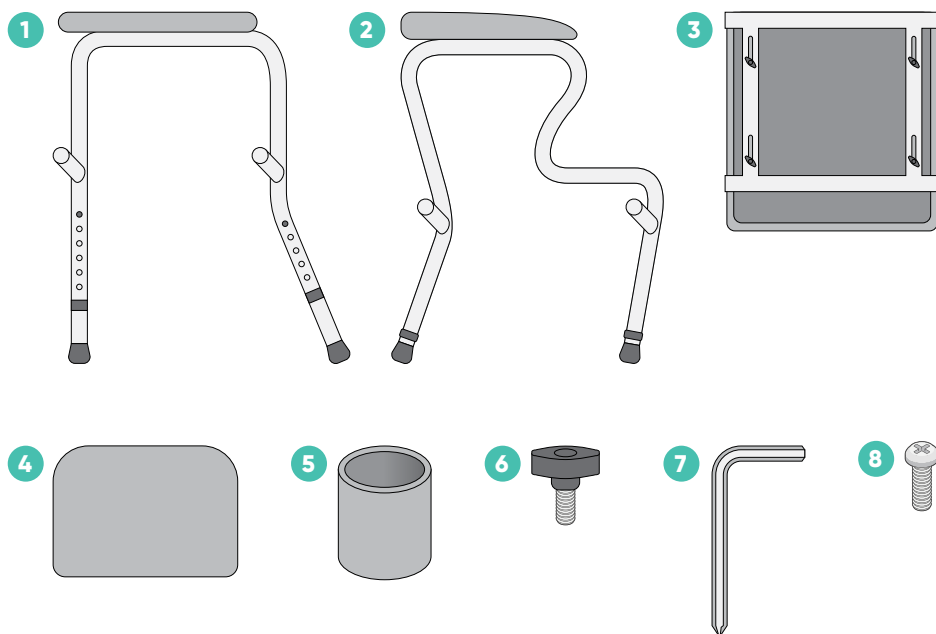

Assembly

01 Setting the Seat Width

Identify which width the seat needs to be assembled to by measuring the user. If the wider seat setting is required, place the spacers onto the below lugs on side frames:

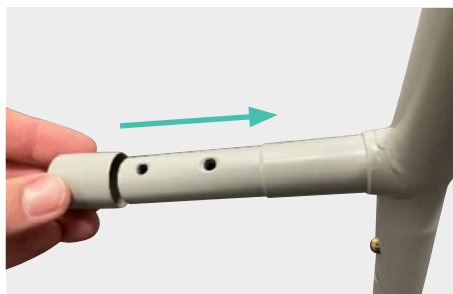

02 Attaching the Seat Frame

Insert side frame lugs into the seat frame until they reach the end.




Assembly

03 Attaching the Padded Backrest

Place the padded backrest against the backrest mounting plates to one side and align the holes (A). Loosely screw the handwheel knobs into the backrest holes (B). Repeat for the other side.

NOTE: Adjust (push or squeeze) the side frames to allow the holes to align. DO NOT TIGHTEN HANDWHEEL KNOBS.

NOTE: the image below shows the backrest being fitted at the wider setting. Holes will align on the outside mounting tab holes if using the narrower setting.

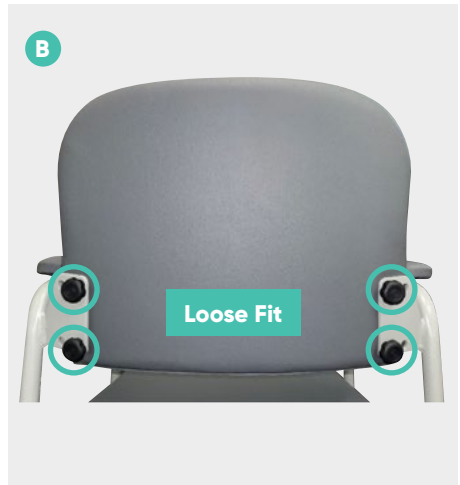

Assembly

04 Securing the Padded Backrest

Secure the seat using the four handwheel knobs (A).

Finally, tighten the backrest handwheel knobs (B).



05 Checking the Push Pins

Ensure push pins on leg adjustment feet are adequately poking out of the tube:



06 Ensure Adequate Fixing

Ensure all fixings are adequately tightened before use.

The chair is now ready to use.

Maximum User Weight

Maximum user weight is 200kg. Ensure the maximum user weight is not exceeded. Please note that maximum weight capacities will be reduced if the user applies their weight unevenly or forcefully.

Heavy or uneven weight distribution may cause stress to the frame and result in damage to the utility chair.

Cleaning & Maintenance

Regular inspection of the desk utility chair for damage to the frame, seat, arms or backrest parts and/or any screws and nuts that may have come loose is recommended.

Use low strength detergent and warm water with a soft textured cloth. Do not use caustic soda, oven cleaner or other high strength cleaning agents.

This product requires no regular servicing aside from cleaning as described above. Please dispose of thoughtfully according to local government regulations.

Warranty Against Defects

Product	Warranty Period
KCare TUBALCO R & R	
All Products (Inc. all electrics)	12 Months
Upholstery	12 months
Shower Chairs, Stools & Toilet Seat Raisers (Zinc Treated products)	36 months structural 5 Years anti-corrosion (external tube only)
Shower Chairs, Stools and Toilet Seat Raisers (Aluminium products)	36 months structural 5 Years anti-corrosion (external tube only)
Shower Chairs, Stools and Toilet Seat Raisers (Stainless steel products)	36 months structural 5 years anti-corrosion
KERRY OXFORD EME EASYVET	
General Equipment	36 months
KERRY	
Manual Handling Equipment	24 months
Slings	12 months
AIRCOMFORT	
Frame Compact Lift Chair	36 months
Frame Day Beds	36 months
Foam	36 months
Air Bags and Gas Struts	12 months
Castors	12 months
Fabric	12 months
Electric motor	24 months
Battery	12 months
MY COMFORT	
MC100 & MC310 Frames	10 years
MC100 & MC310 Electrics	2 years

Tested Against Standards

This product has been tested against and complies with international standard AS/NZS 4688.2:2000 and AS/NZS 4688.3:2000

1. K Care Healthcare Equipment provides consumers with the following warranty in relation to this Product, in addition to complying with the requirements of any relevant legislation, including the Competition and Consumer Act 2010 (Cth) in Australia and the Consumer Guarantees Act 1993 in New Zealand (the Acts), except where a New Zealand consumer acquires the Product for the purposes of a business.
2. In this warranty, we have used the following definitions:
 - (a) Form means a warranty claim form issued by K Care Healthcare Equipment in respect of Products.
 - (b) K Care Healthcare Equipment or our means K Care Healthcare Equipment,
 - (c) Products means the goods manufactured by K Care Healthcare Equipment (including products manufactured by its contract manufacturers) set out in the table above;
 - (d) Material means a material or component used by K Care Healthcare Equipment in the manufacture of the Products;
 - (e) Retailer means the authorised dealer of Products from whom the Product was purchased;
 - (f) Warranty Period means the warranty period set out in the table above, commencing from the date of purchase of the relevant Product. If the Product or part (as the case may be) is repaired or replaced, there will be no extension to the original warranty period; and
 - (g) Workmanship means the handling, assembly and manufacturing processes performed by or on behalf of K Care Healthcare Equipment in order to manufacture the Products.
3. K Care Healthcare Equipment warrants that for the duration of the Warranty Period, all Products will be free of faults arising from defects in Workmanship or Materials, on the terms and conditions set out in this warranty.
4. K Care Healthcare Equipment undertakes that if during the Warranty Period any Product, or any part of a Product, has failed to operate correctly due to faulty Workmanship or defective Material, it will repair or replace the Product or part (as the case may be) free of charge (excluding freight fees if applicable) provided that the following procedure is met:
 - (a) The consumer must contact the Retailer or K Care Healthcare Equipment upon becoming aware of any defect to a Product. The consumer will then be provided with a Form, which must be completed by the consumer and returned to K Care Healthcare Equipment together with satisfactory proof of purchase.
 - (b) K Care Healthcare Equipment will review the Form to determine whether there is a defect, and if so K Care Healthcare Equipment agrees to (at its option) repair, replace or supply equivalent goods, or pay the cost of any of those remedies to the consumer.
 - (c) If K Care Healthcare Equipment requests the return of the applicable Product or part, the consumer will be responsible for the collection and freight costs of returning that Product or part to K Care Healthcare Equipment.
 - (d) The consumer acknowledges that it is also responsible for the freight costs to deliver any new Product or replacement part to it.
5. To the extent that the supply of Products is a supply to a consumer within the meaning of the applicable Act, the consumer is entitled to consumer guarantees which K Care Healthcare Equipment does not exclude, restrict or modify. In all other respects, and to the extent permitted by law, K Care Healthcare Equipment:
 - (a) Limits its liability for any non-excludable condition or warranty to rectifying any defect at its option, as set out in paragraph 4 (b); and
 - (b) Excludes liability (whether express or implied) of any nature whatsoever for any consequential loss, damage or injury arising as a result of any fault in the Products.
6. This warranty does not extend to damage to Products which occurs during transit or transportation, or which is caused by any abuse, accident or improper installation, connection, use, adjustment or repair or use of goods otherwise than in accordance with instructions issued by K Care Healthcare Equipment.
7. The warranty on Products is waived if any addition or attachment to the Products do not have K Care Healthcare Equipment's approval or are not sold as K Care Healthcare Equipment products. The Products are designed to perform specific tasks under established test loads and unauthorised attachments may produce stresses for which the design is not appropriate.
8. The following applies to consumers who purchased a relevant Product in Australia: Our goods come with guarantees that cannot be excluded under the Australian Consumer Law. You are entitled to a replacement or refund for a major failure and compensation for any other reasonably foreseeable loss or damage. You are also entitled to have the goods repaired or replaced if the goods fail to be of acceptable quality and the failure does not amount to a major failure.

Notes

[illegible]

Improving your life everyday™

At Your Service

Sales, administration
and general enquiries

National 1300 783 783

Your Customer Service Team
customerservice@kcare.com.au

kcare.com.au

All information is considered to be true and correct at the date of publication. Changes in circumstances after the time of publication may impact on the accuracy of the information.

