

 **aspire** OVERBED TABLE  
SETUP GUIDE



# 1. OVERBED COMPONENTS

This Overbed Table comes with the following parts:

Description	Quantity
Table Top with bracket mounted on bottom surface	x 1
Upright Post with quick adjustment lever and fixing knob	x 1
C-Shape Base Frame with 4 x black castor wheels (2 x locking)	x 1
14mm Bolts (Long) with washers and spacers (FOR BASE)	x 2
12mm Bolts (Short) with washers (FOR TABLE TOP)	x 2

Thank you for choosing this Aspire Overbed Table.

This product has been designed to offer supreme stability, safety and adjustability.

It is important that you read this setup guide carefully before using your Aspire Overbed Table.

12 Month Warranty

**LAMINATE TABLE TOP** accommodates food trays, personal items, and provides a stable writing surface

**FIXING KNOB** to tighten upright in place

**QUICK ADJUSTMENT LEVER** makes height adjustment easy for users and carers

**FLAT PACK CARTON DESIGN** is utilised to improve shipping / storage efficiency and reduce potential for damage

**TWO LOCKING CASTORS** allow the table to move easily and be positioned securely



Laminate AC  
White Top  
BEA006700



Laminate  
Beech Top  
BEA006701



Thermoform Recessed  
Top (Cream)  
BEA028000



Split Tilt Top  
(Cream)  
BEA028100

## 2. ASSEMBLY INSTRUCTIONS

Open the carton carefully. DO NOT USE a box cutter or sharp implement as this may damage the product. For OH&S, lay the carton on its side with the top of the carton open and slide out the Overbed Table contents.


1. Insert spacers, 14mm bolts and washers into the mounting plate on the C-shape base (from below).



2. Locate 14mm bolts into the base of the upright post and tighten. Ensure the upright post is straight on the base plate whilst tightening. Loctite or other fixing agents may be used if required.



3. Insert upright post table support bar into table top bracket and attach using 12mm short bolts and small washers.

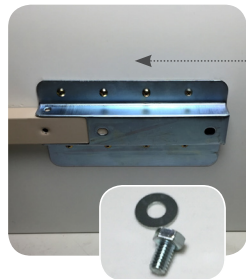

4. Remove the transportation screw from the upright and your over bed table is ready for use.



## 3. ADJUSTING YOUR TABLE

To Adjust your Over Bed Table:

1. Loosen the fixing knob on the upright post.



2. With one hand pull and release the quick adjustment lever and move the table up or down until the desired height is reached. Minimal force should be required.



3. Re-tighten fixing knob to lock in place.



**Imported by:**

AIDACARE P/L | 4 Noonan Rd, Ingleburn NSW 2565

T - 1300 133 120 | W - [www.aidacare.com.au](http://www.aidacare.com.au)

MADE IN CHINA