


AGED CARE & HOME CARE BEDS



“...strong local businesses that have a broad product portfolio...”

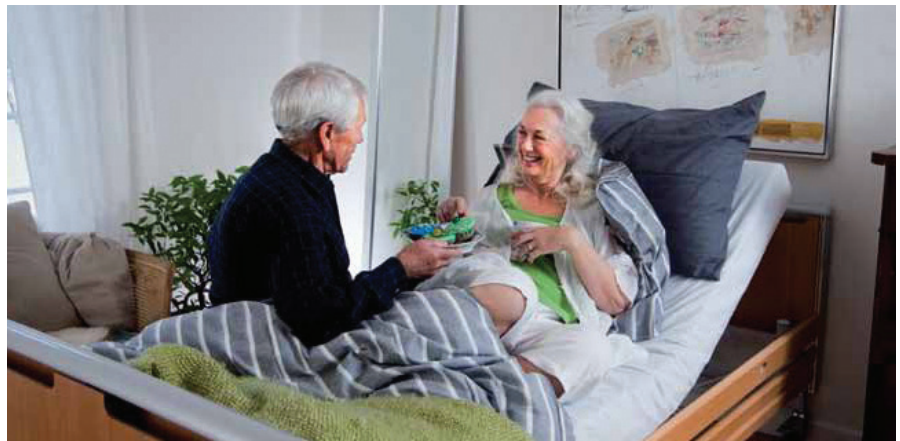

125 years of experience

Invacare is the world's leading manufacturer of home healthcare and mobility equipment. We proudly supply Australian homes, nursing homes, clinics, and hospitals with innovative solutions to promote recovery and active lifestyles. With the largest distribution network in the industry, servicing more than 80 countries, we are truly a global leader. But how did we get to where we are today?

The roots of Invacare Corporation stretch back to 1885, when Ohio tricycle company Fay Manufacturing (later Worthington Manufacturing) began production of bicycle-wheeled mobility chairs for persons with disabilities. These were of course the early versions of modern-day wheelchairs and would go on to form the core of the business. In the 1950s, the company became 'Mobilaid', and was still mainly a wheelchair manufacturer when it was bought by Boston Capital Corporation (BCC) in 1971 and merged with Invalex Company, to form Invacare Corporation.

Today, Invacare is the world's top manufacturer of rehabilitative products, respiratory products, home care beds, soft goods & supplies, institutional beds & equipment, and of course, wheelchairs.

Invacare Australia





Invacare has grown both organically and through the acquisition of Australian organisations, to create strong local businesses that have a broad product portfolio to support institutional rehabilitation and home health care.

Invacare Australia's team of professionals work closely with dealers and clinicians nationwide to service customers with their rehabilitation, home healthcare and mobility equipment needs.



“...quality products, focused on rebuilding active lifestyles...”
for disabled and aged consumers...



					
	1600	CS7	CS7 Wide	Etude® Plus	Octave®
Care Environment					
Home Care				✓	✓
Aged Care	✓	✓	✓		✓
Hospital Care	✓	✓	✓		✓
Sleep surface functions					
Electric backrest	✓	✓	✓	✓	✓
Electrical knee break	✓	✓	✓	✓	✓
Electric leg rest					✓
Electric height adjustment	✓	✓	✓	✓	✓
Tilt/Trendelenberg		✓ Option	✓ Option	✓	
Braking System					
Toe brake				✓	
Central braking	✓	✓	✓		✓
Specifications					
Sleeping surface width	800 mm	900 mm	1054 mm	900 mm	1190 mm
Width overall	920 mm	900 mm	1054 mm	1020 mm	1190 mm
Sleeping surface length	1900 mm	2032 mm	2032 mm	2000 mm	2000 mm
Length overall	2139 mm	2140 mm	2140 mm	2190 mm	2210 mm
Under bed clearance	130 mm	356 mm (centered for 2/3)	356 mm (centered for 2/3)	190 mm (minimum)	150 mm
Bed height range	200 to 670 mm	178 to 762 mm	178 to 762 mm	330 to 730 mm or 400 to 800 mm	400 to 800 mm
Caster size	125 mm	50 mm	50 mm	100 mm	125 mm
Safe Working Load (SWL)	200 kg	226 kg	226 kg	180 kg	385 kg
Maximum Patient weight (35 kg for acc)	165 kg	190 kg	190 kg	145 kg	350 kg
Electronics	LINAK IPX4	LINAK IPX4	LINAK IPX4	LINAK IPX6	LINAK IP66
Compatible Bed ends					
Specific Bed ends				✓	✓
Queen Anne, Arch, Colonial, Standard	✓				
Amherst - Solar Oak		✓	✓		
Drake - Biltmore Cherry		✓	✓		
Custom Bed ends and surrounds	✓	✓	✓		
Compatible Side Rails					
Steel Side rails	✓	✓		✓ Options	
Wooden Side rails				✓	✓
Assist Rail		✓	✓		
Assist Bar		✓	✓		
Support Handle		✓	✓		✓
Self Help Pole	✓	✓	✓	✓	✓
Accessories					
Bumper bars	✓	✓	✓		
IV pole	✓			✓	✓
Mattress extension	✓	✓	✓	✓	✓
Available	AU	ANZ	ANZ	ANZ	ANZ


Guide to Understanding Features

Invacare with its aim to making it simple for you to use our catalogue and products has provided a very simple key that gives you information at your fingertips.


Beds




Mattress Base /
Overall Length

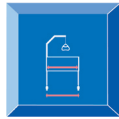


Mattress Base /
Overall Width

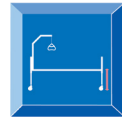

Mattress Base Height
Adjustment

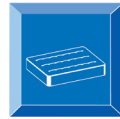


Recommended
Mattress Size

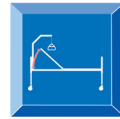


Back Rest
Adjustment




Leg Rest Lift

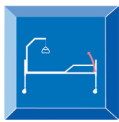

Thigh Angle

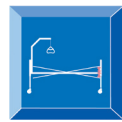

Tilt Option

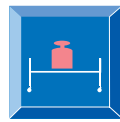

Safe Working
Load*

* Safe Working Load - 35kg = Suggested User Weight, in accordance with EN1970, based on an allowance of 35kg for accessories, including the mattress.

INVACARE®

1600




From traditions of Murphy Furniture and Hencicare Hospital Equipment comes the 1600 aged care bed, designed specifically for the Australian market. Tried and trusted, it remains the first choice for many carers today.

Options

- Range of head and foot board designs. Custom timber surrounds are also available.
- BZA-H151220RG Bed extension (25cm)
- BSH-3093 Metal horizontal safety side rails
- BZA-H2133 IV pole
- BZA-H2116 Self help pole
- BZA-B161205RG Wall protection
- Buffer bar


Central locking brake system

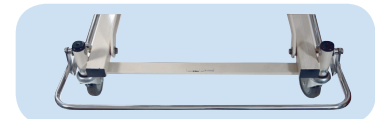

Model PKG-BAS-1601-C125RG

Features	Benefits
1600 mattress base height can be lowered to 20cm	Reduces risk of injury and enhances patient independence
4 section mattress support	Knee break to suit individual patient requirements
Choice of head and foot board design and finish	Offers variety and flexibility
Sheet metal base	Perfect for use in aged care facilities
Patient lifter clearance of 13cm	Suitable for use with most patient lifters, for ease of use with new or existing equipment

Technical Specifications

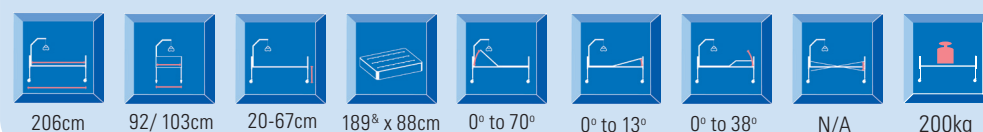

⁸Standard configuration only, when fitting a mattress for use with a 1600, use a longer model (189 x 88cm) or add a bolster.

INVACARE® CS5



Contoured slat deck design improves infection control, and reduces mattress wear and tear.

Intelligent design ensures that the CS5 is perfectly vertical when raising or lowering, with no tracking or travelling taking place.

The CS5 bed was developed for our customers requiring mobility at any height. The trademark CS Series features tool-less assembly and is easy-to-clean.

CS Accessories



Side rails (IHCSRLPA)



Self Help Pole (IHCSTRPHD)



Assist bar (IHCSBAAS)
Also available: 10cm bed extension (IHCSBEK)



Lowers to 23cm from the ground



Key lock system secures the head and foot board in place

Model IHCS5-CL

Features	Benefits
Tool-less assembly	Fast connect system allows for faster set-up and easier alterations
Centrally locked rear castors	Floor lock system secures the bed to the floor at all height ranges but allows for mobility on demand
Standard 4 section with knee break	For maximum comfort and care
Weight capacity of 226kg	High weight capacity ensures that the CS5 is suitable for a wide variety of residents
Range of accessories available	Accessories include metal side rails, assist bar and lifting pole - all of which are interchangeable between the CS Series of beds



Amherst style head and foot boards in Kensington Maple (IHAM)

Technical Specifications

212/219cm	90cm	23/76cm	197 x 88cm	0° to 70°	N/A	0° to 30°
						226kg

INVACARE® CS7



Contoured slat deck design improves infection control, and reduces mattress wear and tear.

Intelligent design ensures that the CS7 does not track or travel during height adjustments, helping minimise wall damage.

An Integrated Attendant Control Panel (ACP) is available and also provides ease for the carer.

The Invacare CS7 Bed is a truly versatile bed. Lowering to a height of 17cm's off the ground, it is a hi-lo bed combined with floor level benefits, making this a completely unique aged care bed, suitable for a variety of residents and their specific needs.



Lowers to a height of 17cm

Options

- Standard bed ends available in Amherst (Solar Oak) IHAM/WCP-TREND or Drake (Biltmore Cherry) IHDR/WCP-TREND. Custom design also available upon request.
- Custom timber surrounds are also available.
- Accessories include Metal side rails IHCSRLPA, Self help pole IHCSSTRPHD, Assist Bar IHCSBAAS
- Pan deck is available.
- Mattress Extension IHCSBEK
- Range of mattresses available to suit individual care needs

Model IHCS7/WS

Features	Benefits
Auto-Contour that simultaneously raises both head and knees	This not only assists in comfort and positioning, but also helps reduce painful skin shearing
Height adjustable between 17-75cm	Lowering to 17cm from the floor, the CS7 is low enough to the floor to help minimise injuries sustained when falling out of bed
Tough and durable hand control pendant	A rugged, water-resistant pendant with fewer recessed areas is easy to clean and very intuitive for the resident to use
Centrally locked rear castors	Floor lock system secures the bed to the floor at all height ranges but allows for mobility on demand
Tool-less assembly and Fast connect system	This allows for faster and easier bed set-up and easier alterations to be made




Drake style head and foot boards in Biltmore Cherry (IHDR/WCP-TREND)



Amherst style head and foot boards in Solar Oak (IHAM/WCP-TREND)

Technical Specifications

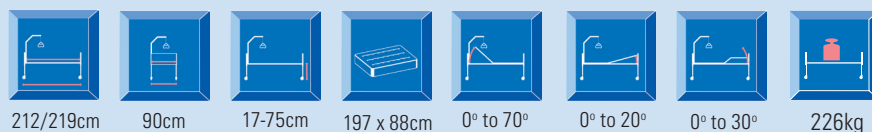

212/219cm

90cm

17-75cm

197 x 88cm

0° to 70°

0° to 20°

0° to 30°

226kg

INVACARE® CS7 Wide Deck

The CS7 Wide Deck bed incorporates the same innovative features as the CS7 with a 42" wide mattress support



42" wide deck offers additional comfort for a variety of residents

Tool less assembly and fast connect system allow for faster set-up and easier alterations

An Integrated Attendant Control Panel (ACP) is available and also provides ease for the carer

Options

- Accessories include 1/2 Rail Assist IHCSRLAS-42, Self help pole IHCSTRPHD, Assist Bar IHCSBAAS-42.
- Pan deck is available.
- Standard bed ends available in Amherst (Solar Oak) IHAM42/WCP-TREND or Drake (Biltmore Cherry) IHDR42/WCP-TREND. Custom design also available upon request.
- Custom timber surrounds are also available.
- Mattress Extension IHCSBEK
- Range of mattresses to suit individual care needs.



Slat deck as standard for improved infection control



Centrally locked low profile castors minimise trip hazards

Model No IHACS42

Features	Benefits
Wider deck/mattress support	At 42" in comparison to the standard 36", the wider deck is prepared for larger residents while also increasing comfort
Height adjustable between 17-75cm	Lowering to 17cm from the floor, the CS7 is low enough to the floor to help minimise injuries sustained when falling out of bed
Safe Working Load of 226kg	With a SWL of 226kg, the CS7 Wide Deck is a bed highly suitable for larger residents
Auto-Contour that simultaneously raises both head and knees	This not only assists in comfort and positioning, but also helps reduce painful skin shearing
Tough and durable hand control pendant	A rugged, water-resistant pendant with fewer recessed areas is easy to clean and very intuitive for the resident to use




Drake style head and foot boards in Biltmore Cherry (IHDR42/WCP-TREND)



Amherst style head and foot boards in Solar Oak (IHAM42/WCP-TREND)

Technical Specifications

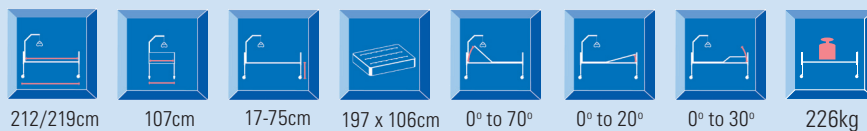

212/219cm

107cm

17-75cm

197 x 106cm

0° to 70°

0° to 20°


0° to 30°

226kg

* May be increased with the optional bed extension/s

INVACARE® Etude Plus

The Etude Plus has been designed to complement furniture in the home whilst offering a multitude of optional accessories to cater for a broad range of users.

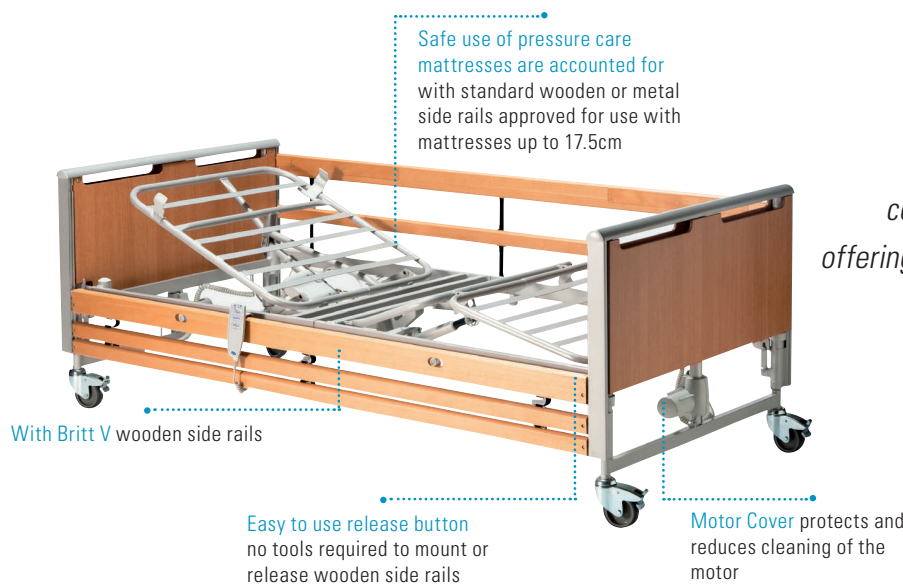

Options

- 1528930-0154 Scala 2 Basic 3/4 steel side rails
- 1538458-0101 Britt V Full length wooden side rails
- 1522786-0125 Aluminium Line Beech Full length side rails
- 1530817-0154 Scala Decubi 2 Foldable 3/4 steel side rails
- 1542093-0154 Self help pole with adjustable handle
- 1543779-0154 Etude Plus transport kit
- 1423994-0154 20cm mattress support extension
- 1529665-0154 Turnable support handle



With metal side rails




In Transport Trolley


ETUDE_4D1544211_PO

Features	Benefits
European design and style designed with smooth soft lines with choice and flexibility of upgrades	All the functions of a nursing home bed which fits in beautifully with other furniture in the home environment
Easy to use hand control operates electrical height adjustment from 40cm to 80cm or 33cm to 73cm through the standard feature of built in double brackets	Provides user with autonomy and assists with caregiver's workload
Auto-regression	Previously not available in the Etude series, auto-regression helps the backrest swing up and back as the bed profiles, helping to reduce shearing
Unique low position of the Britt V wooden side rails (option, as shown above)	A low position of just 3cm above the mattress support level on the accessed side ensures minimum squeezing for the user when leaving or entering the bed
Tool-free construction	Transport, assembly & disassembly by one person
Meets the new European safety standards	Etude Plus bed and accessories are compliant with the updated IEC 60601-2-52

Technical Specifications

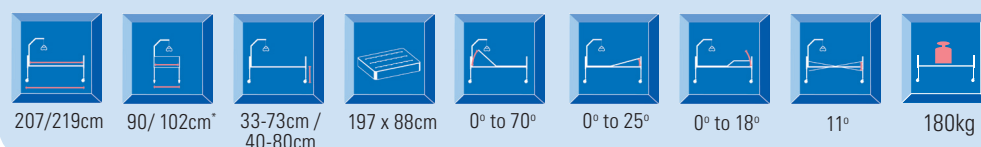


INVACARE® Octave

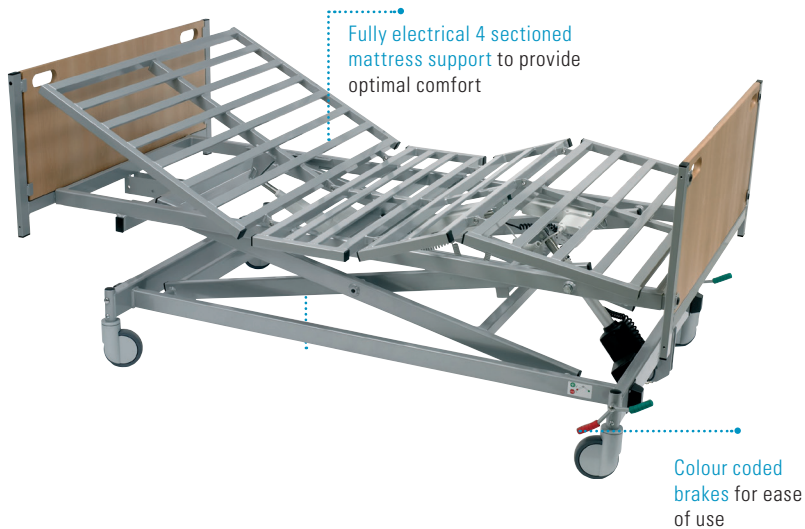

The Octave bed has been designed with the bariatric client in mind with a safe working load of 385kg. The frameless design ensures optimum patient care, further backed by its state of the art electrical profiling system.

Options

- 1446353-0152 Support handle
 - 14463500-0101 Wooden adapter - kit for fitting wooden side rails
 - 1427875-0101 Wooden side rails
 - 50.57700.MO Lifting pole#
 - 1478730-0152 Transport Trolley#
 - 10 cm extension* 1446498-0152
 - 1451259 Decapo mattress#
- # Included in Full Kit *Up to 2 can be fitted

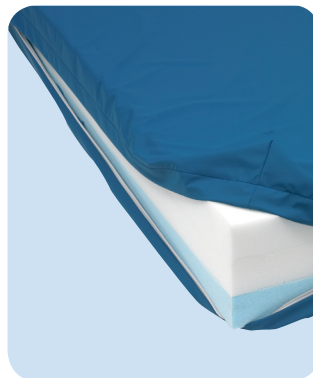

Decapo mattress




With Decapo mattress, lifting pole and side rails

Model OCTAVE_C2AB1_M1_1

Features	Benefits
Safe Working Load of 385kg and User Weight 350kg	For those who require additional support
Simple and easy to use hand control unit operates electrical height adjustment from 40cm to 80cm	Provides user with autonomy and assists with carers' workload
Easily detachable bed ends	The carer has the best possible access to the user
"Frameless" mattress support with bed frame that extends to the edges of the bed	The mattress sits on the edge of the bed frame enabling easier access for the carer and better support for the user. Easy transfer of patient in and out of bed
Fully electrical 4-sectioned mattress support as standard	Providing optimal comfort for the user
Transport Trolley as standard	Delivery, transportation and storage of the Octave are made as convenient as possible

Technical Specifications

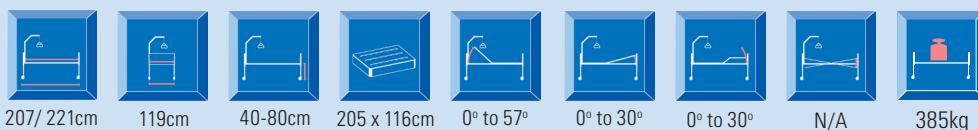

INVACARE® Furniture

The Invacare range of furniture is designed to match any aged care surrounding. From over bed tables, to custom timber bedside lockers, Invacare offers a wide range of options and configurations to suit your needs.



C- BASE OVER BED TABLE

Technical Data

Standard height: 72cm - 110cm

Powder coated steel frame

'C'-shaped base

Durable melamine surface with PVC edges

Twin wheel castors

Standard colours: Wiluna White, Cherry, Select Beech and Champagne

STANDARD COLOURS

Wiluna White, Cherry, Select Beech and Champagne

WARDROBES

The Invacare wardrobe line includes a wide range of door and shelf configurations, designed to suit the storage requirements of any resident or facility. The custom wardrobe range is also available in a wide variety of colours.



Colonial Bedside

Cambridge Bedside

BEDSIDE CABINETS

Invacares® bedside cabinets come in a full range of drawer and cupboard combinations to suit individual storage requirements. They are available in a wide selection of colours and handle styles that will compliment the room setting. A choice of two or four twin wheel castors is provided for mobility and cleaning. Towel rail and locking top drawer are available options.



Laura Mitred Arm Chair



Manuela Full Back Chair

DINING TABLES AND DAY CHAIRS

Our range of Australian made custom dining tables and day chairs will help complete any aged care setting and is available in a variety of styles, configurations and finishes. They are also manufactured from top grade material, ensuring that quality is to the highest standards.



“ When you combine the power of **innovation** and **care**, the result is **Invacare.** ”



Invacare Australia Pty Ltd
 ABN 45 074 676 378
 PO Box 5002, 1 Lenton Place,
 North Rocks, NSW 2151,
 Australia.
 Freephone: 1800 460 460
 Fax: 1800 814 367
 Email: orders@invacare.com.au
 Web: www.invacare.com.au

Invacare New Zealand
 PO Box 62-124, 4 Westfield Place,
 Mt. Wellington, Auckland,
 New Zealand.
 Freephone 0800 468 222
 Freefax 0800 807 788
 Email sales@invacare.co.nz
 Web: www.invacare.co.nz

Invacare Asia Ltd.
 PO Box 5002, 1 Lenton Place,
 North Rocks, NSW 2151,
 Australia.
 Phone: +61 2 8839 5330
 Fax: +61 2 8839 5331
 Email: exportsales@invacare.com.au
 Web: www.invacare.com

For more pre-sales information about this wheelchair, please consult the applicable website.
 For information about typical user description, intended use, adjustments and transportation, please refer to the product user guide.
 PB7196/1214



Yes, you can.